

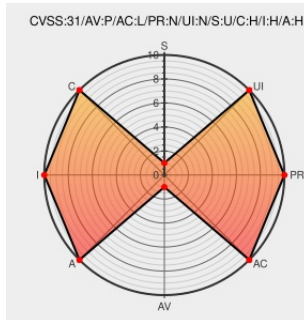


CVE-2019-13546

Published on: 10/25/2019 12:00:00 AM UTC

Last Modified on: 03/23/2021 11:27:45 PM UTC

CVE-2019-13546

[Source: Mitre](#)[Source: NIST](#)[CVE.ORG](#)[Print: PDF](#) 

Certain versions of [Intellispace Perinatal](#) from [Philips](#) contain the following vulnerability:

In IntelliSpace Perinatal, Versions K and prior, a vulnerability within the IntelliSpace Perinatal application environment could enable an unauthorized attacker with physical access to a locked application screen, or an authorized remote desktop session host application user to break-out from the containment of the application and access

unauthorized resources from the Windows operating system as the limited-access Windows user. Due to potential Windows vulnerabilities, it may be possible for additional attack methods to be used to escalate privileges on the operating system.

CVE-2019-13546 has been assigned by ics-cert@hq.dhs.gov to track the vulnerability - currently rated as **MEDIUM** severity.

CVSS3 Score: **6.8 - MEDIUM**

Attack Vector	Attack Complexity	Privileges Required	User Interaction
PHYSICAL	LOW	NONE	NONE
Scope	Confidentiality Impact	Integrity Impact	Availability Impact
UNCHANGED	HIGH	HIGH	HIGH

CVSS2 Score: **7.2 - HIGH**

Access Vector	Access Complexity	Authentication
LOCAL	LOW	NONE
Confidentiality Impact	Integrity Impact	Availability Impact
COMPLETE	COMPLETE	COMPLETE

CVE References

Description	Tags	Link
-------------	------	------

By selecting these links, you may be leaving CVEreport webspace. We have provided these links to other websites because they may have information that would be of interest to you. No inferences should be drawn on account of other sites being referenced, or not, from this page. There may be other websites that are more appropriate for your purpose. CVEreport does not necessarily endorse the views expressed, or concur with the facts presented on these sites. Further, CVEreport does not endorse any commercial products that may be mentioned on these sites. Please address comments about any linked pages to comment@cve.report.

There are currently no QIDs associated with this CVE

Known Affected Configurations (CPE V2.3)

Type	Vendor	Product	Version	Update	Edition	Language
Application	Philips	Intellispace Perinatal	All	All	All	All
cpe:2.3:a:philips:intellispace_perinatal:*:*:*:*:*:*:						

No vendor comments have been submitted for this CVE

[← Previous ID](#)

[Next ID →](#)

© [CVE.report](#) 2023   |

Use of this information constitutes acceptance for use in an AS IS condition. There are NO warranties, implied or otherwise, with regard to this information or its use. Any use of this information is at the user's risk. It is the responsibility of user to evaluate the accuracy, completeness or usefulness of any information, opinion, advice or other content. EACH USER WILL BE SOLELY RESPONSIBLE FOR ANY consequences of his or her direct or indirect use of this web site. ALL WARRANTIES OF ANY KIND ARE EXPRESSLY DISCLAIMED. This site will NOT BE LIABLE FOR ANY DIRECT, INDIRECT or any other kind of loss.

CVE, CWE, and OVAL are registered trademarks of [The MITRE Corporation](#) and the authoritative source of CVE content is [MITRE's CVE web site](#). This site includes MITRE data granted under the following [license](#).

CVE.report and Source URL Uptime Status [status.cve.report](#)